

Cleaning and Disinfecting Guide

Disinfectant List for DR Systems

- The following disinfectants can be used to disinfect appearance & cover of digital X-ray systems
: GC85A, GM85, GC80, GC70

Name	Manufacture	Active Ingredient	CBC*	EPA*
Accel PREvention Wipes	Diversey Canada, Inc.	Hydrogen peroxide		
ALI-FLEX RTU	Lalema	Sodium hypochlorite		
ANIOSURF Premium NPC	ANIOS	Quaternary Ammonium(N-Alkyl)		
ANIOS社 surfa'safe	ANIOS	Quaternary Ammonium(N-Alkyl)		
Cleanisept Wipes MAXI	Dr.Schumacher GmbH	Quaternary Ammonium(N-Alkyl)		
CLOROX Bleach Germicidal	The Chlorox Company	Sodium hypochlorite	○	○
Dettol	Reckitt Benckiser	Quaternary Ammonium(N-Alkyl)		
Ecolab Actichlor plus	Ecolab	Troclosene sodium		
ED Wipes	MH healthcare	Quaternary Ammonium(N-Alkyl)		
Enzysept PLUS Wipes	Wool C&Tech	Quaternary Ammonium(N-Alkyl)		
Incidin OxyWipes	Ecolab	Hydrogen peroxide		
Metricide-OPA plus	Metrex	Ortho-phthalaldehyde		
Perasafe	LANXESS Deutschland GmbH	Citric acid		
Protex Ultra	Parker Laboratory Inc.	Quaternary Ammonium(N-Alkyl)		
SaFecide + (Sensitive wipes 200)	Liebecos	Quaternary Ammonium(N-Alkyl)		
Sani Cloth Active wipes	PDI	Quaternary Ammonium(N-Alkyl)		
Sani Cloth Universal wipes	PDI	Quaternary Ammonium(N-Alkyl)		
Sani-cloth bleach	PDI	Sodium hypochlorite	○	○
Sani-Cloth HB	PDI	Quaternary Ammonium(N-Alkyl)		
Sodium hypochlorite	Yuhan Chlorox	Hypochlorous acid		
Sono Ultrasound Wipes	Advanced Ultrasound Solutions Inc.	Quaternary Ammonium(N-Alkyl)	○	
Sporox II	DSHealthcare Inc.	Hydrogen peroxide		
Tristel duo	Tristel Solutions Limited	Chlorine dioxide		
Tristel Sporicidal Wipe	Tristel Solutions Limited	Chlorine dioxide		
T-spray II	Pharmaceutical Innovations, Inc.	Quaternary Ammonium(N-Alkyl)		
Virkon S	LANXESS	Pentapotassium, Succinic acid		
Wip Anios Excel	Laboratoires Anios	Quaternary Ammonium(N-Alkyl)		

※ Recommended by CDC (Centers for Disease Control and Prevention) for USA Market.
CBC : Center for Biocide Chemistries.
EPA : Environment Protection Agency


Cautions & Warnings for disinfecting DR system

- Use only cleaning agents and disinfectants included in the list, especially for marked areas.
- Rub the product carefully with soft cloth when using disinfectants.



- Be careful not to let liquids or other substances get inside
: Do not spray directly on the product, use a cloth or tissue.

- ※ Wear a face mask and gloves when cleaning or disinfecting the product.
- ※ Do not clean any electrical connectors, it can cause an electric shock or fire.

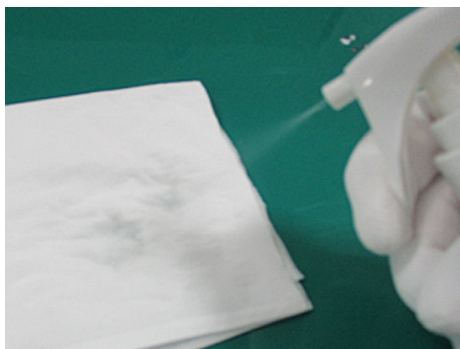

Clean the surface in accordance with the disinfectant manufacturer's instructions (Temp, wiping/contact time etc.)

Disinfectant List for Detector All

· The following disinfectants can be used to disinfect detectors

: S-Detector™ (S4343-W, S4335-W, S3025-W), AccE Detector (S4343-AW, S4335-AW),
AccE Standard Detector (S4343-AWM, S4335-AWM), AccE Glass-Free Detector (F4335-AW)

Name	Manufacture	Active Ingredient	CBC*	EPA*
Accel PREvention Wipes	Diversey Canada, Inc.	Hydrogen peroxide		
ALI-FLEX RTU	Lalema	Sodium hypochlorite		
ANIOSURF Premium NPC	ANIOS	Quaternary Ammonium(N-Alkyl)		
ANIOS社 surfa'safe	ANIOS	Quaternary Ammonium(N-Alkyl)		
Anioxyde 1000	Laboratoires Anios	Peracetic acid		
Azomax	Vernacare	5-chloro-2-methyl-2h-isothiazol-3-one		
Cleanisept Wipes	Dr.Schumacher GmbH	Quaternary Ammonium(N-Alkyl)		
Cleanisept wipe forte	Dr.Schumacher GmbH	Quaternary Ammonium(N-Alkyl)		
Cleanisept Wipes MAXI	Dr.Schumacher GmbH	Quaternary Ammonium(N-Alkyl)		
Clinell Sporidical wipes	GAMA Healthcare Ltd	Sodium Percarbonate		
Clinell Universal Wipes	GAMA Healthcare Ltd	Quaternary Ammonium(N-Alkyl)		
CLOROX Bleach Germicidal	The Chlorox Company	Sodium hypochlorite	○	○
CN'F	CARROLLASIA Co.,LTD	Alcohol		
Dettol	Reckitt Benckiser	Quaternary Ammonium(N-Alkyl)		
Distel Disinfectant Wipes	Tristel	Polymeric Biguanide Hydrochloride		
Ecolab Actichlor plus	Ecolab	Troclosene sodium		
ED Wipes	MH healthcare	Quaternary Ammonium(N-Alkyl)		
Enzysept PLUS Wipes	Wooil C&Tech	Quaternary Ammonium(N-Alkyl)		
Incidin OxyWipes	Ecolab	Hydrogen peroxide		
Lysoform	Lysoform	Formaldehyde		
Metricide-OPA plus	Metrex	Ortho-phthalaldehyde		
Mikrobac Tissues	BODE Chemie GmbH (HARTMANN)	Quaternary Ammonium(N-Alkyl)		
Perasafe	LANXESS Deutschland GmbH	Citric acid		
Protex Ultra	Parker Laboratory Inc.	Quaternary Ammonium(N-Alkyl)		
SaFecide + (Medical wipes 200)	Liebecos	Alcohol		
SaFecide + (Sensitive wipes 200)	Liebecos	Quaternary Ammonium(N-Alkyl)		
Sani-cloth bleach	PDI	Sodium hypochlorite	○	○
Sani-Cloth HB	PDI	Quaternary Ammonium(N-Alkyl)		
Sani-Cloth Super	PDI	Quaternary Ammonium(N-Alkyl)		
Sodium hypochlorite	Yuhan Chlorox	Hypochlorous acid		
Sono Ultrasound Wipes	Advanced Ultrasound Solutions Inc.	Quaternary Ammonium(N-Alkyl)	○	
Sporox II	DSHealthcare Inc.	Hydrogen peroxide		
Tristel duo	Tristel Solutions Limited	Chlorine dioxide		
Tristel Sporidical Wipe	Tristel Solutions Limited	Chlorine dioxide		
Trophon Companion wipe	Nanosonics Limited	Quaternary Ammonium(N-Alkyl)		
T-spray II	Pharmaceutical Innovations, Inc.	Quaternary Ammonium(N-Alkyl)		
Virkon S	LANXESS	Pentapotassium, Succinic acid		
Wip Anios Excel	Laboratoires Anios	Quaternary Ammonium(N-Alkyl)		

※ Recommended by CDC (Centers for Disease Control and Prevention) for USA Market.


CBC : Center for Biocide Chemistries.

EPA : Environment Protection Agency

Cleaning Guide for S-Detector™ & AccE Detector


Cleaning

- Always turn off the equipment completely before cleaning.
- Do not insert liquids into the system.
- When cleaning the exterior of the system, use a soft cleaning cloth dampened with tepid Water and soap.

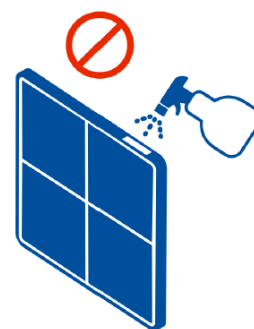


Liquid intrusion

- Do not spill liquid or chemicals onto the instrument or, in cases where the patient is injured, do not allow it to become wet with blood or other body fluids, as doing so may result in fire or electric shock.
- In such situation, protect the instrument with disposable covering as necessary.

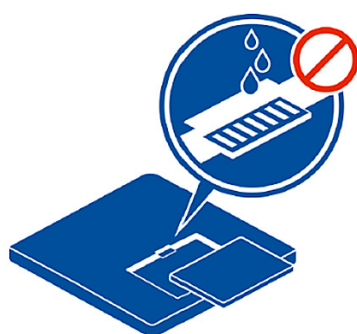


S-Detector™ only

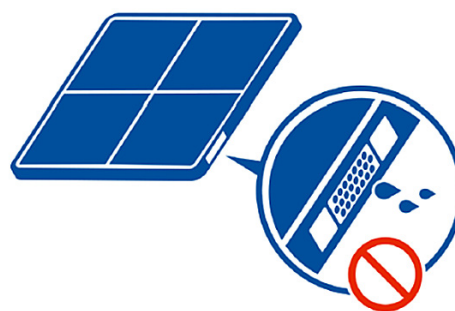

Cleaning Guide for AccE Detector

- Recommended using a commercial plastic brush when cleaning non-slip pads.
- If need, non-slip pad can be exchanged or detached.

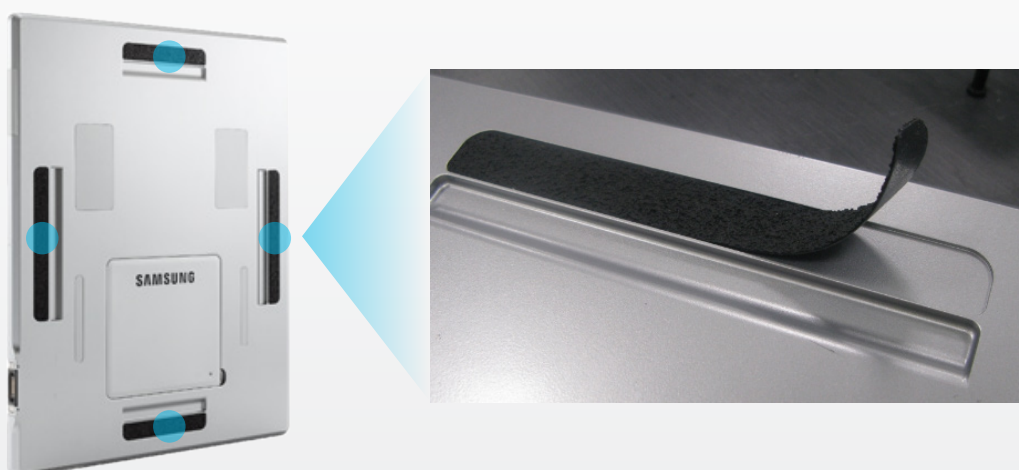

Usage Guide for Detector Storage

- Detector has to be equipped properly.
- If you equip the detector in the storage rear for charging or moving with too much force or diagonally, pogo pin on the main cable connector may get damaged.
- There shouldn't be any moisture on the detector connector before equipping detector.
- If you equip the detector on the storage rear before completely removing moisture after cleansing, the moisture may get on the main cable connector, resulting in the damage of pogo pin.

